

CONTACT ARENA TO COORDINATE ORDERS

☎ 602-357-3288

THE arena[®] TOUR MASK

THE FIRST AND ONLY AMERICAN
MADE, WASHABLE, REVERSIBLE,
N95 FILTER COMPATIBLE FACE MASK

OUTER MATERIAL SPECS

OUTSIDE SHELL COMPOSITION:

100% Polyester

COLOR:

Black, White, Navy,
Dark Navy, Grey, Tan

CORPORATE BRANDING:

Full color printing on front or
back panels (*Now Available*)

* REVERSIBLE WHITE

WEIGHT:

250 gsm

THICKNESS:

3 MM

OPTIONS:

Top & Bottom Tie
Back or Top Tie
Back w/ Bottom
Elastic Strap

* MADE IN THE USA



*CUSTOM PRINT SAMPLES
AVAILABLE FOR PURCHASE.

onemillion
ENTREPRENEURS

PROUD SUPPORTER

arena® TOUR MASK

The Tour Mask is constructed using a 2-piece athletic-grade composition, with a reversible design. The 4-way stretch fabric supports the body's natural cooling process by wicking away moisture evenly throughout the construction of the fabric. The outside polyester shell maintains the shape and provides an additional layer of filtration for large particles.

- We use 3M Proprietary N95 microfiber filter media, designed for easy breathing.
- N95 Certified 3M Filter Media is cut to fit edge-to-edge of each mask size.
- Includes Aluminum Nose Clip for a more formed seal.

PM 2.5 CARBON FILTER SIZE: 12CM * 8 CM / 4.72IN * 3.15 IN
3M N95 FILTER MEDIA SIZE: 24.8CM * 15.6CM / 9.75IN * 6.15IN

WASH & REUSE UP TO 50X

The athletic-grade fabrics and elastics allow for heavy machine washing and industrial sanitizing. We recommend washing the masks prior to first use. Having two masks per operator is optimal for daily "wash-one, wear-one" applications.

The interior shell contains a pocket that can hold the optional 3M Proprietary Filter Media or a PM2.5 Carbon Filter.



- * REVERSIBLE WHITE
- * FULL COLOR IMAGE

INTERIOR MATERIAL SPECS

INTERIOR SHELL COMPOSITION:

82% Poly / 18% Spandex

COLOR:
White / Full Color

WEIGHT:
240 gsm

THICKNESS:
0.35MM

ADDITIONAL INFO:

- Ties or Elastic Straps
- 4 way stretch fabric
- High-performance fibers that support the body's natural cooling process by wicking away moisture evenly throughout the construction of the fabric
- Peach skin finish for a soft hand feel throughout the inner shell.

3M PROPRIETARY FILTER MEDIA:

The mask is designed to hold a compatible 3M proprietary filter media full-mask insert or a PM2.5 disposable carbon filter. Inserts can be replaced through the side seams of the mask and center into place between the vertical darts on the nose bridge and chin. The inner shell is the pocket for the disposable inserts.



PRICING TIERS


Mask Exterior Shell with/without Insert
(Available to ship now):

QUANTITY	PER UNIT	N95 FILTER	CUSTOM PRINT
1-499	\$12.99	\$1.00	\$3.00
500-999	\$11.99	\$1.00	\$3.00
1,000-99,999	\$10.99	\$1.00	\$2.00
100,000-249,999	\$8.99	\$1.00	\$1.50
250,000-499,999	\$7.99	\$1.00	\$1.00
500,000-1MM+	\$7.49	\$1.00	\$0.75

arena®


Tour Mask

This is a quick usage guide. Please read the most recent instructions & FAQs online: tourmask.com

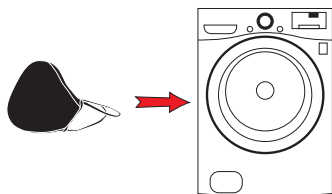

WASH BEFORE USE

Wash before first use. We recommend laundering the masks daily after use. The materials we've used can withstand 190F temperatures. DO NOT wash the filter. Always remove the filter before washing.

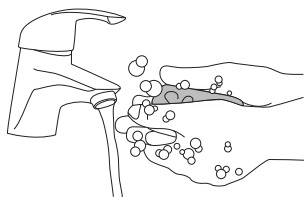


WASH YOUR HANDS


Always wash your hands with soap and warm water for at least 20-seconds - before AND after you handle the mask.

INSERT FILTER (optional)


The Tour Mask is compatible with standard N95 / PM2.5 filters. Filters are disposable, and should be changed frequently.

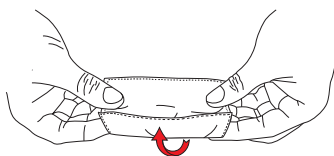


1) To Insert the disposable filter, roll the top and bottom edges of the filter media into a cylinder.

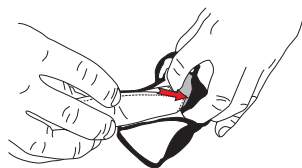

2) Carefully begin to insert the filter through the side opening, on the inside of the mask.

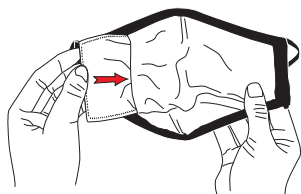

3) Continue until the filter is centered between the top nose and bottom chin darts.

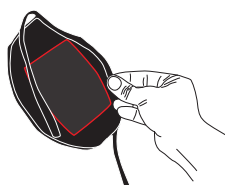

4) Check to ensure the filter is correctly positioned in the mask.


FITTING MASK



1) Place the mask so that it covers your nose and mouth. We recommend that the mesh side faces outwards.



2) Pinch the shapeable aluminum nose strip to help form a better seal under the eyes and around your nose bridge.

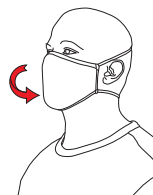

REVERSING THE MASK

Reversing the mask (so the white side faces outwards) can help with breathability - for extended daily wear. We DO NOT recommend this option for use in public places. Go online for how to reverse the mask.

REMOVING MASK




1) Wash your hands - then carefully pull the straps over your head. Use the straps to pull the mask away from your face.




2) Never touch the mask as you risk contaminating yourself - keep your hands off of your face and use common sense.

REMOVE FILTER

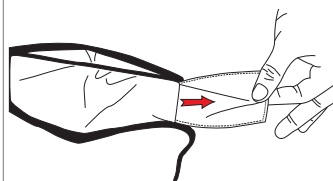

1) Change the filter as frequently as possible, ideally daily. The filter can last for up to 7-days depending on the manufacturer.

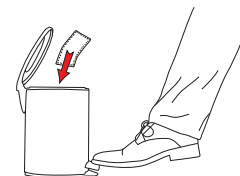

2) Always discard the filter in a safe place immediately after removing it from the mask - do not leave the filter lying around.